

Positioning with a difference

Belts and harnesses from Posigo and Soloc




Posigo FLO2 Anterior Trunk Supports

Breathe freely

Traditional anterior trunk supports mainly focus on positioning the torso, often compromising between support and comfort.

The Posigo FLO2 anterior shoulder harness is a welcome exception to this: the innovative FLO2 system provides maximum support across the shoulders and upper torso without restricting the user's breathing.

Dynamic straps placed below the sternal support pad allow the harness to move in concert with the diaphragm, enabling unrestricted breathing.

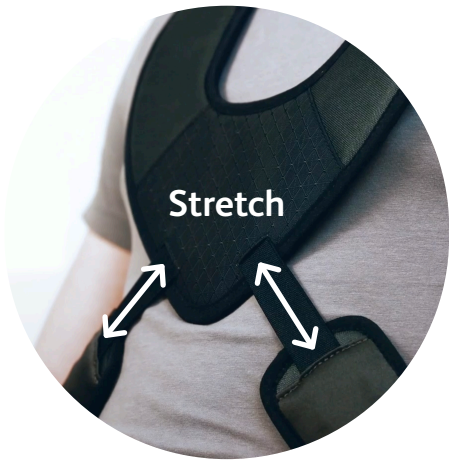

FLO2



See more here



Posigo Pelvic Positioning Belts

Modular 4-1 solution

The Posigo selection of pelvic positioning belts offers a comfortable and secure option for positioning the seated person safely in their wheelchair. They provide a modular 4 in 1 solution with unique features and benefits:



Adjustability

Each pad on the Posigo pelvic positioning belt has a rotatable knob to adjust the belt length and tension, allowing adaptability for different clothing on different days. The knobs are positioned so that they can be adjusted by either the occupant or the carer.



Interchangeable buckle

The standard push button buckle on the Posigo pelvic positioning belt can be interchanged with other buckles, such as a security buckle that needs a device such as a pen to push through the reduced diameter hole.



Attachable secondary straps

Posigo pelvic positioning belts have optional secondary straps which clip onto each pad to create a 3-point or 4-point belt. This makes the belt easy to adapt appropriately for the individual needs around different pelvic positioning challenges.

See more here



Mounting guides

The Posigo pelvic positioning belts are supplied as 2-point belts with full mounting kits. The mounting straps are provided with numbered markings to facilitate specifying the mounted strap lengths on each side.

Meet safety standards

Posigo pelvic positioning belts are tested to AS EN ISO wheelchair seating strength and flammability standards.



Soloc Pelvic Positioning Belts

Three Belts. One Magnetic Buckle. Total Support.

For those with limited hand dexterity

1

Soloc Connect

The standard support belt for everyday positioning and security.

2

Soloc Uno

Designed for one-handed use – easy to fasten and release with a single hand.

3

Soloc Freedom

Ideal for users with limited or no hand movement – operates with minimal effort.



The buckle

Effortless Magnetic Buckle

All Soloc belts feature an easy-to-use magnetic buckle – perfect for users with limited hand function, limb difference, or one-handed use.

Secure by Design

Magnets guide the connection, while a mechanical lock keeps the belt safely closed, even under pressure.

Simple to Release

Push the button or pull the tab – the buckle opens smoothly, even with full body weight against it.

See video here



Soloc Connect

The Connect is a standard pelvic positioning belt with the Soloc buckle attached. It allows users with reduced hand function to fasten and release their own belt. This reduces reliance on carers and increases independence.

The belt is adjustable. You can add padding to aid with postural positioning. Padding is available in two different lengths (600mm or 800mm).



Soloc Uno

A first in pelvic positioning belts – the Uno is designed for those with the use of only one hand/arm. One side of the buckle is held in place by a FlexiFix arm – positioned ready for connection. Great for stroke survivors or those with hemiplegia.

Simply pull the stud side across to meet the buckle – the magnets make attachment easy. Then pull the belt tight using the large D-ring. You can fit the belt for left or right handed use. Padding is available in two different lengths (600mm or 700mm).



Soloc Freedom

The Freedom is an innovative belt that has been designed for a unique user group. It can be fastened by pushing down on the FlexiFix arms. This enables those with no hand function or even those with no hands at all to use it independently!

The belt stays in an open position for transfers in and out of the chair – no more sitting on the belt and being unable to release it!



Posigo Dynamic Foot Supports

Unobtrusive functionality

The Posigo Dynamic Foot Supports use elasticated straps to provide stability during sudden movements, helping the user stay well positioned while maintaining freedom of movement. The supports also enable dynamic positioning of the foot on the wheelchair foot support, with a magnetic closure for easy on-and-off use. As a result, users enjoy greater autonomy.

Placement

Attachment of a foot support at the ankle encourages plantar reflexes which stabilise the foot on the foot support.

Ergonomic design

The Posigo ankle support has been designed to fit around the ankle area for maximum support and comfort.

Dynamic strapping

The securing straps are all elasticated to allow a controlled amount of movement for the foot.

Magnetic buckle

The length of the strap around the ankle can be adjusted using hook and loop attachments. The fastening buckle is a magnetic closure, siting and attaching the buckle very easily.

Foot support attachment

The straps for attaching the Posigo ankle support to the foot support have numbered gradations to aid in achieving the required lengths. Each pair of straps comes with a mounting kit for attachment to the foot supports.

Safety

The Posigo ankle supports are tested to the AS EN ISO wheelchair seating standards.

Choice of sizes

The Posigo ankle support comes in five sizes from XS to XL.

